

Dondolino makes the active inclusion of the child in his surrounding environment easier when used in the standard quality of life







the experience of manual activities guided by vision improves manipulative activities and the eye-hanc coordination



anti-gravity muscle practice, which is the recovery of a partial anti-gravity setting of posture



facilitation of the circulatory functions thanks to the orthostatic position that helps a more effective bloodstream



facilitation of peristalsis thanks to the orthostatic position that helps







news!

also size 1 with trunk support

height and rotation adjustment with the new "safe lock" system

technical features



tilting frame for children who need to improve their righting reflex. Continuous adjustment, without snaps (by gas spring)



no edge frame, in order to enable trunk rotation and side-arm motion



"TOOL FREE" adjustments by knobs, without using wrenches



graduated bars to keep the postural setting chosen by the therapist



birch wood table. It is removable and adjustable in height tilt and depth



birch wood foot board. It gives a safe guide to the child It is adjustable in height



twinned solid pivoting wheels, 75 mm. diam., rear wheels with brakes. Thanks to the wheels the aid can be moved from one room to another having the child on it


adjustments


multi-adjustable trunk and pelvic supports, they are padded, wide and with breathable upholstery


ergonomic knee-pad supports

heel rests

accessories

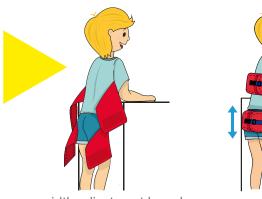

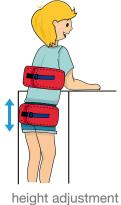


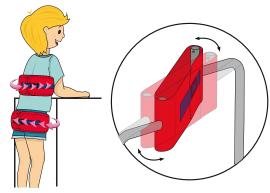


rotation adjustment, to adapt to the posture of the child






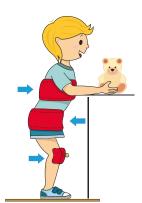


independent back and forth adjustment. should the child have enough trunk control the trunk support can be removed

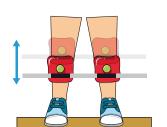


height adjustment

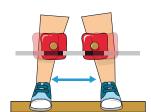


width adjustment

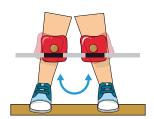

rotation adjustment

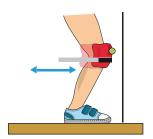

depth adjustment

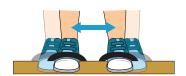

width adjustment



depth adjustment




944 circumference reduction parts both for trunk and pelvic support. Thickness: 3 cm. each side



865 adjustable headrest height, tilt and depth adjustable



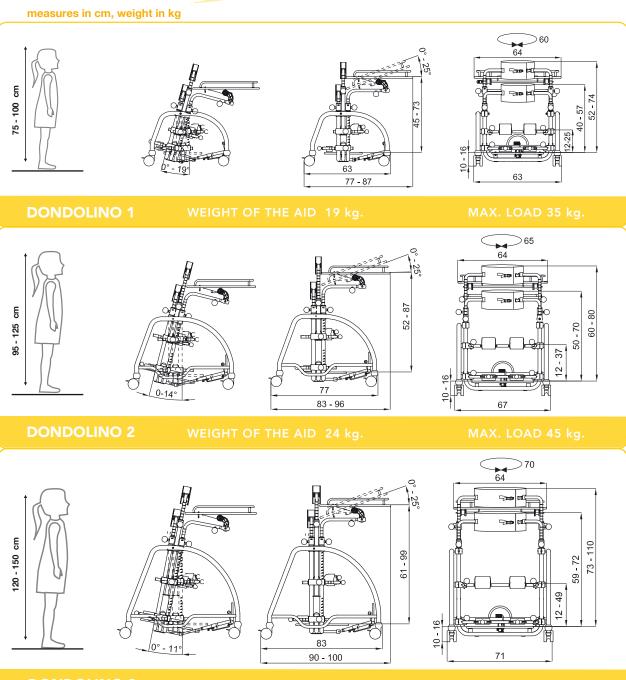

hypen-separated numbers indicate the minimum and maximum range of adjustment

PRODUCT MARKED **(**

Authorized Dealer

The copyright even if partial is not allowed without the authorization of the manufacturer. The manufacturer reserves the right to make changes without any notice



Via A. Da Sangallo,1 - 06034 FOLIGNO (PG) ITALY TEL. +39 0742 22927-FAX +39 0742 22637 info@ormesa.com-www.ormesa.com



