Ò.



Convaid Flyer[®]

Lightweight, Modular and Portable Tilt-in-Space Wheelchair

Convaid

Convaid Flyer[®]

The Convaid Flyer is a lightweight, modular and portable wheelchair developed with youth-oriented aesthetics and growth in mind. A width-expandable chair that a child can grow with, the Flyer has a uniquely designed tilt-in-space base that also accepts more complex adaptive seating. The Flyer comes standard with -5° to +40° of tilt, and 30° of back angle adjustment, offers optional recline and a large selection of options and accessories to achieve ideal positioning and support.

The modular design allows for easy removal of the seating module from the mobility base for simple lifting and compact transport. The Flyer seating module has the ability to grow in width from 12" to 16" and seat depth can grow 4", providing therapists and caregivers the convenience and ease of use they desire without sacrificing the seating and positioning needs the user requires.



Convaid

Lightweight, Modular and Portable

- Lightweight, Modular and Portable Removal of the lightweight seating module (16.8 lbs.) from the mobility base assists in easy lifting, transport and storage.
- Growth-adaptable Seating module can grow 12-16" in width, and 4" in depth; total seat depth range is 11-21".
- Tilt and Recline Comes standard with 45° of tilt, 80-110° of back angle adjustment and an optional 80-170° of recline.



Emmy with mom and her Convaid Flyer®.

SEATING MODULE OPTIONS



Back Canes



Solid Back Seating System





7.5" x 2" Front Caster, 16" Rear Tire

(17" Front Seat to Floor Height)

Seating Module

TRANSIT



Meets WC19 Safety Standards for Safe Transportation


TIRE OPTIONS





6" x 2" Front Caster, 12" Rear Tire (15" Front Seat to Floor Height)

FRAME COLOR OPTIONS

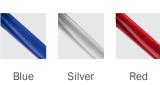

UPHOLSTERY ACCENT COLOR OPTIONS

Mobility Base



Please visit Convaid-R82.com to view all accessories offered with the Flyer®.

Please scan OR code to visit website and view additional product information



CONVAID FLYER® SPECIFICATIONS

	Seat Size	12	13	14	15	16
	Seat width	12"	13"	14"	15"	16"
/\	Minimum seat width with hip guides	9"	10"	11"	12"	13"
<u> </u>	Seat depth options	Total available seat depth range is 11"-21"; customer to select range setting based on back cane position; each seat depth has 4" of growth.				
<u>VI</u>	Seat back height options	3 back cane height options to achieve 19", 22" and 25" With back pans, range is 13"-25"				
<u>م</u>	Seat to footplate standard 90° fixed	Range is 6"-11.5"				
<u>6</u> *	Seat to footplate elevating legrest	9-14" 5.25-11.25"*	9-14" 5.25-11.25"*	11-16" 7.25-13.25"*	11-16" 7.25-13.25"*	13-18" 9.25-15.75"*
<u>921</u>	Seat to floor	2 wheel options provide 15" and 17" STF				
P.	Seat to back angle	80-110° (standard)				
F.	Back recline	80-170° (option)				
Ś	Tilt range	-5° (forward) to 40° (rear)				
6	Overall height 15"STF/17"STF	34-40"36-42"Height varies based on selected back cane height and selected seat to floor height				
<u>Þ</u> ,	Overall length	25.75" or 28" Varies based on selected tires				
	Overall width	21"	22"	23"	24"	25"
	Shoulder strap height	16-24" range with back pans				
5	Push handle adjustments	±9" after selecting back cane size				
ibs i	Weight capacity (regular)	Back Cane Position 1: 85 lbs. Back Cane Positions 2-4: 170 lbs.				
	Weight capacity (transit)	Back Cane Position 1: 85 lbs. Back Cane Positions 2-4: 170 lbs.				
Â	Weight of chair (varies based on tires)	12" tires: 30.7 lbs. (mobility base = 13.9 lbs., seating system with leg rest = 16.8 lbs.) 16" tires: 34 lbs. (mobility base = 17.2 lbs., seating system with leg rest = 16.8 lbs.)				
* with flipped footplates						

Everyone should be able to live a free and independent life and pursue their dreams regardless of any physical circumstances.

Convaid Products, LLC 2830 California Street Torrance CA 90503 Tel: 1-844-US-Mobility (844-876-6245) Fax: 310-618-2166 www.Convaid-R82.com CB2001 0218

