

## The New Convaid Trekker



## Content


- 1 The New Convaid Trekker
- 2 Convenient, Lightweight, Compact
- 3 Adaptive Seating
- 4 Standard Features
- 5 Product Overview
- 6 Trekker Accessories
- 8 Specifications



### The New Convaid Trekker®

Convaid teamed up with R82 to bring you our newest generation tilt-in-space wheelchair, covering a wide range of users, from early intervention up to 110 lbs. The new Trekker has an improved docking interface, making it easy for you to slide the seat and lock it when you need to attach, detach or reverse the seat. Fold your Trekker frame easily into a compact and lightweight package with or without the seat and use only one hand to unfold.

Enjoy the highly versatile seating options offered with new Trekker: Depending on your child's need during the day, you can use the Trekker seating module with the R82 High Low indoor frame or use the R82 Stingray size 1 seat with the Trekker frame.

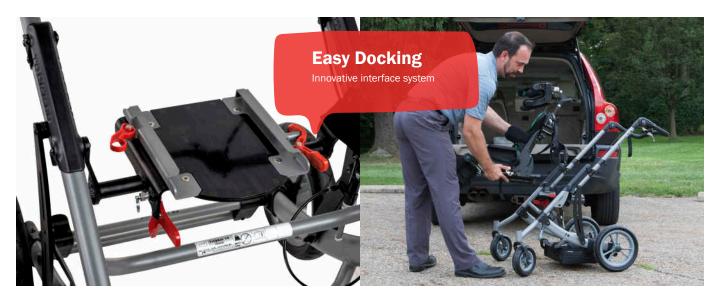







## **Convenient, Lightweight and Compact**

The Trekker is very lightweight and its seat and frame fold down to a low height making it easily portable. The Trekker can fold with or without the seat and only requires one hand to unfold.





## **Adaptive Seating**

Research suggests that adaptive seating systems may be able to improve activity and participation at home among children with severe cerebral palsy (Angusupailsal 2015). The Trekker can be used to establish an ideal seating position for children and youngsters with disabilities to play and engage in activities that help develop their skills.



#### A new way for children to be involved without compromising comfort and interaction

For most children and youngsters, a supportive seating position can have a positive impact on muscle tone and facilitate the child's ability to move their head with their upper extremities. A range of different features and functions have been developed to create the best possible seating position for each individual child.

#### Now Compatible with R82 Stingray Seat and High-Low Indoor Frame

The Convaid Trekker has a NEW docking interface that allows for use of an extended seating range including the R82 Stingray Size 1 seating module, as well as the attachment of the Trekker seating module to the R82 High Low Indoor Frame. The Trekker has always been and remains a light weight, compact folding Tilt in Space Wheelchair. It is a best-in-class product and can be customized to each individual's physical needs. The Trekker is made to last and designed to provide the optimal comfort, safety and mobility for the child.



Trekker Seating Module on Trekker Frame



R82 Stingray Seating Module on Trekker Frame



Trekker Seating Module on High-Low Indoor Frame



R82 High-Low Frame



### **Standard Features**

- Adjustable Tilt in Space, -10°-45°
- Seat to back angle, 80° 110° recline
- 180° Reversible seating
- Foot activated wheel locks
- Quick release rear wheels
- Separate frame & seating module
- 90° Fixed front rigging
- 90° Fixed footplate
- One piece height adjustable push handle



**RECLINE AND TILT** 

The clinical benefits of 45° of posterior tilt are promotion of proper pressure relief, weight distribution, and reduction in the risk of skin breakdown.



MULTIPLE GROWTH OPTIONS

Seat depth, height, width adjustable. Firm laterals for thoracic support. Hip guides are width, height and depth adjustable for pelvic support.



The Trekker can be positioned carriage or stroller style



**MULTI-ADJUSTABLE** 

Height adjustable seat pan to accommodate up to 26 inches of seat to back height.



**REVERSE SEATING** 

Reversible so you can have your child face forward or face you.









\*Positioning options & accessories sold separately




## **Trekker Accessories**









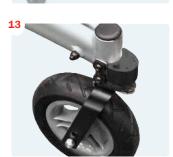






б





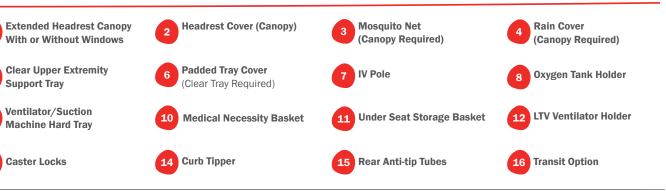

Seating | Trekker

Convaid | R82 | a part of Etac



Seating | Trekker

7



## Compatible with R82 Stingray Seat Size 1

- Easy tilt-in-space and back recline adjustments
- Postural support and stability
- 7 fabric colors made of 100% recycled material



# R82 Stingray Seating

Specifications

Seat Size	1
Seat width	6.25 - 11.75"
Seat depth	7-11.75"
Back height	16-21.5"
Seat width, folded	16.5"
Seat length, folded	22.5"
Seat height, folded	11"
Lower leg length	5.5-17.5"
Leg rest angling	-30° -> 30°
Back recline	170°
Weight, seat	12 lb
Max. user weight/ load, seat	88 lb
Max. user weight/load, transportation	88 lb



Convaid Trekker Seat			Convaid Trekker Frame		
Seat Size	1 (12")	2 (14")	Frame Size	1 (14")	
Seat width	11.75"	13.75"	Folded length	39"	
Seat depth range	7.5" - 11.5"	9.5" - 16.5"		26"	
Back height	13" - 17"	17" - 26"	Folded width		
Hip width between hip guides	6.5"-11.5"	6"-14"	Folded height	13.5"	
Trunk width	8.5"-12"	10.5"-14"	Push handle adjustments	29" - 43"	
Trunk width between	6.5"-12"	8.5"-14"	Weight of frame	25 lbs	
swing-away laterals	4" - 12	4" - 12"	Overall height	37.5"	
Seat to Footplate (one-piece footplate)	3" - 11"*	3" - 11"*	1	42.5"	
Seat to Footplate (individual footplates)^	4.75" - 10.5" 1" - 5"*	7" - 13" 1" - 6.25"*	Overall length	42.5	
Back angle adjustment	80° - 110°	80° - 110°	Overall width	26"	
Optional recline adjustment	80° - 170°	80° - 170°	Weight capacity regular/ transit	110/110 lbs	
Tilt adjustment	-5° to 45°	-5° to 45°			
Headrest extension	8"	8"			
Shoulder strap height	14.5" - 20.5"	15" - 23"			
Footplate adjustment angle	±15°	±15°			
$\fbox Weight of seating system^{\Delta}$	10 lbs	11.5 lbs			
			*Reduced height when footplates are flipped ∆Weight of chair is without front rigging		